Patient Support Rails and Patient Support Rails Plus

[English](#)

[Deutsch](#)

[Français](#)

[Español](#)

[Dansk](#)

[Suomi](#)

[Italiano](#)

[Svenska](#)

[Norsk](#)

[Nederlands](#)

[Polski](#)

[Türkçe](#)

For Model:

9A600001

9A600002

User Guide

Style U

003-10114-99 Rev. AA9 (10/11/23)

Product Information

Dealer:

Date of Purchase:

Model / Serial Number:

Midmark Authorized Service Company:

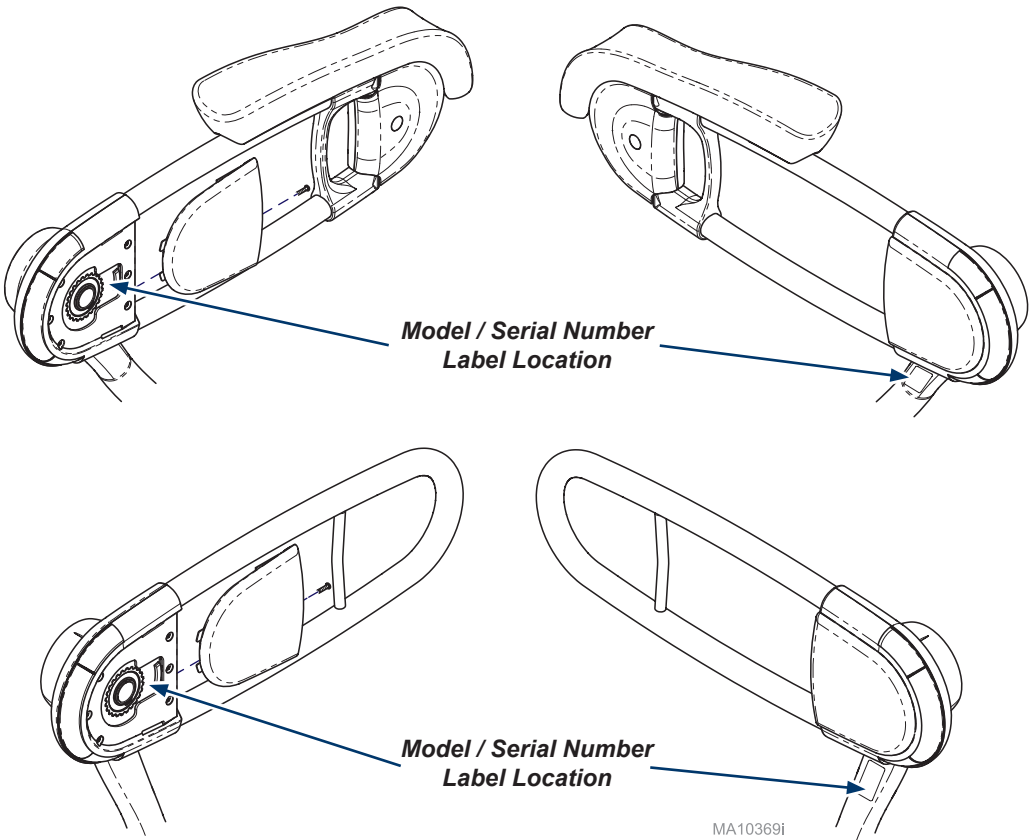
Table of Contents

<u>Product Information</u>	2
<u>Product Registration</u>	4
<u>Transportation / Storage Conditions</u>	4
<u>Operating Conditions</u>	4
<u>Disposal of Equipment</u>	4
<u>Authorized Representatives</u>	4
<u>Important Information</u>	6
<u>Intended Use</u>	7
<u>Safety Instructions</u>	7
 Operation - Patient Support Rail	
<u>Lowest Position</u>	8
<u>15° Patient Position</u>	8
<u>30° Patient Position</u>	9
<u>Access Position</u>	9
 Operation - Patient Support Rail Plus	
<u>Lowest position</u>	10
<u>15° Patient Position</u>	11
<u>30° Patient Position</u>	12
<u>Access Position</u>	13
<u>Ingress Egress Position</u>	13
<u>15° Phlebotomy Articulation</u>	14
<u>30° Phlebotomy Articulation</u>	14
 Dimensions	
<u>Seat Width</u>	15
 Maintenance	
<u>Calling for Service</u>	16
<u>Cleaning</u>	16
<u>Periodic Maintenance</u>	16
 Warranty Information	
<u>Limited Warranty</u>	17

Product Registration

SCAN ME

Transportation / Storage Conditions

Ambient Temperature Range:.....-30°C to +60°C (-22°F to +140°F)

Relative Humidity..... 10% to 90% (non-condensing)

Operating Conditions

Ambient Temperature Range:.....+10°C to +40°C (+50°F to +104°F)

Relative Humidity.....30% to 75% (non-condensing)

Disposal of Equipment

At the end of product life, the chair, accessories, and other consumable goods may have become contaminated due to normal medical use. Consult local codes and ordinances for proper disposal of equipment, accessories and other consumable goods.

Authorized Representatives

EC

REP

WARNING

*In the event of any serious incident in relation to the device,
please contact Midmark and the appropriate competent authority.*

Customers in the EU should direct all questions, incidents and complaints to the Midmark Authorized Representative listed below.

CEpartner4U
Esdoornlaan 13
3951 DB Maarn, The Netherlands
Phone: +31 343 442 524
Fax: +31 343 442 162

Customers in the UK should direct all questions, incidents and complaints to the Midmark UK Responsible Person listed below.

UKCApartner4U Ltd.,
7 CAMPION WAY, BINGHAM, NOTTINGHAM NG13 8TR, UK.
www.UKCApartner4U.co.uk
Phone : +44 7905 384429

Customers in Australia should direct all questions, incidents and complaints to the Midmark Sponsor listed below.

ICONA Pty Ltd.
1A/2A Westall Road
Clayton, Victoria 3168
Phone: 1300 442 662

003-10114-99

English - 4

© Midmark Corporation 2020

TP202 20-42-FO-00014 Rev A1 C2169

Authorized Representatives - continued

Customers in India should direct all questions, incidents and complaints to the Midmark Authorized Representative listed below.

Midmark India
Kalpataru Point, Unit No. 12, 1st Floor Kamani road, Sion (E)
Mumbai, 400022
Tel: +91 22 4915 3000
Fax: +91 22 4915 3100

Customers in Saudi Arabia should direct all questions, incidents and complaints to the Midmark Authorized Representative listed below.

Attieh Medico Ltd
Abdullah M. Al Khalifi Street
Al Nakheel Dist. II
P.O. Box 116105
Jeddah (East)
21391
Saudi Arabia
Tel: + 966 2 286 4707
Fax: + 966 2 286 4744

Customers in the United Arab Emirates should direct all questions, incidents and complaints to the Midmark Authorized Representative listed below.

GULF DRUG LLC.
Al Barsha,
Dubai, UAE
Phone: +971 4 501 4000
Fax: +971 4 501 4100

Customers in Hong Kong should direct all questions, incidents and complaints to the Midmark Authorized Representative listed below.

Associated Medical Supplies Company LTD. 沒有中文版使用說明
Room 1201-1202, 12/F, Fo Tan Industrial Center
26-28 Au Pui Wan Street, Fo Tan, Shatin
New Territories, Hong Kong
Phone: +852 2604 9389
fax: +852 2694 0866

Customers in Qatar should direct all questions, incidents and complaints to the Midmark Authorized Representative listed below.

Tadmur Trading
Midmac Roundabout
Opposite Al Mana Petrol Station
Doha, Qatar
Tel: +974.44337000
Fax: +974.44337100

Customers in Israel should direct all questions, incidents and complaints to the Midmark Authorized Representative listed below.

Medtechnica Ltd.
HaTnufa St 7
Petah Tikva, Israel
Tel: +972 3.925.4040
Fax: +972 3.924.9977

Important Information - Safety Symbols

WARNING
Indicates a potentially hazardous situation which could result in serious injury.

Caution
Indicates a potentially hazardous situation which may result in minor or moderate injury. It may also be used to alert against unsafe practices

Equipment Alert
Indicates a potentially hazardous situation which could result in equipment damage.

Note
Amplifies a procedure, practice or condition.

Symbol Glossary

Pressure Limit

Temperature Limit

Proper shipping orientation

Keep dry

Fragile

Refer to manual

Manufacturer

Catalogue Number

Serial Number

Warning: Crushing of Hands

Medical Device

Intended Use

Patient Support Rail: The patient support rail is an ADA compliant product that is intended to assist patient access to the chair and assist the patient in repositioning while on the chair.

ADA Side Rail with Articulating Arm Pad: The side rail with articulating arm pad is an ADA compliant product that is intended to assist patient access to the chair and assist the patient in repositioning while on the chair. In addition it is intended to properly position a patient's arm for blood pressure checks and phlebotomy procedures.

Safety Instructions

WARNING

No modification of this equipment is allowed.

WARNING

Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

WARNING

Do not sit on patient support rail or armrest pad.

WARNING

Do not press any button during a procedure.

WARNING

Do not use arm pad as a gripping surface during ingress egress.

WARNING

Avoid pinch point under arm pads.

WARNING


Avoid entrapment between rail and back.

Equipment Alert

42lbs is the recommended maximum pressure on the armrest pad.


Operation - Patient Support Rail

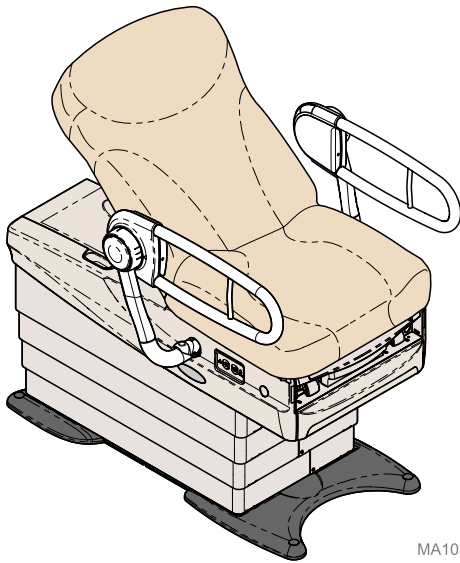MA10356i

Lowest position


Starting position of Support Rail.

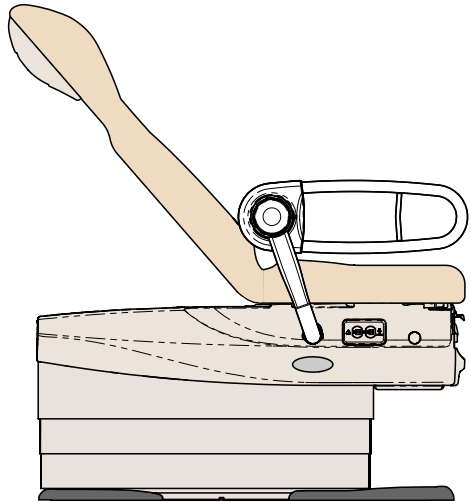MA10357i

Patient Right - View


This position is used for arm positioning with short patients.

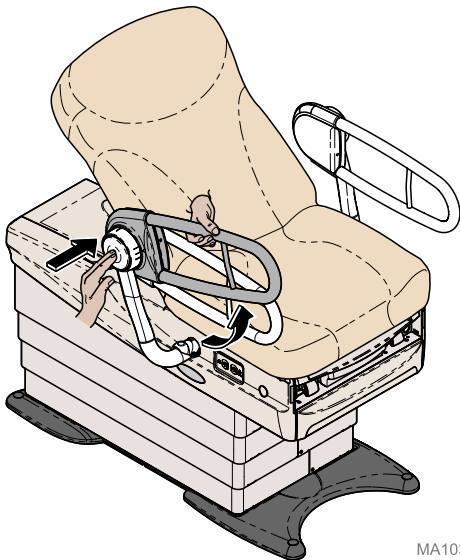MA10358i

Articulation of rail

Hold rail with hand while pushing button with other hand to raise support rail into position.

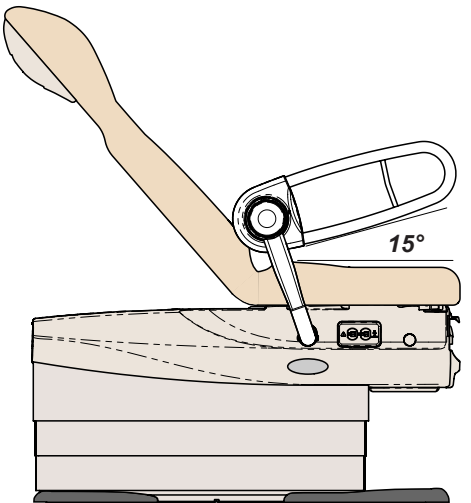
MA10359i

Patient Right - View


The 15° position is used for arm positioning with medium height patients.

Operation - Patient Support Rail

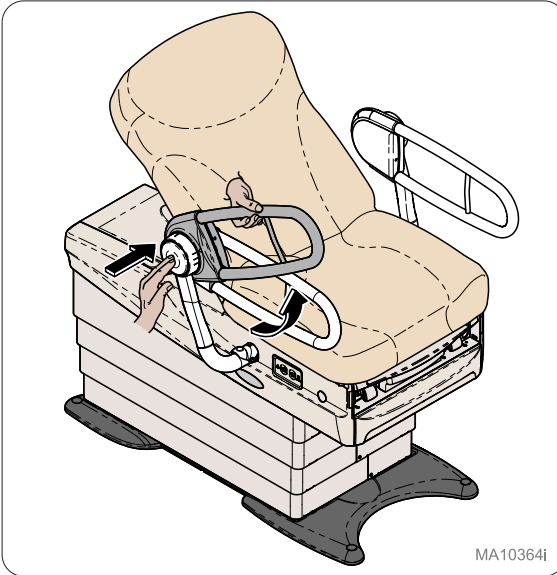
Articulation of rail

Hold rail with hand while pushing button with other hand to raise support rail into position.

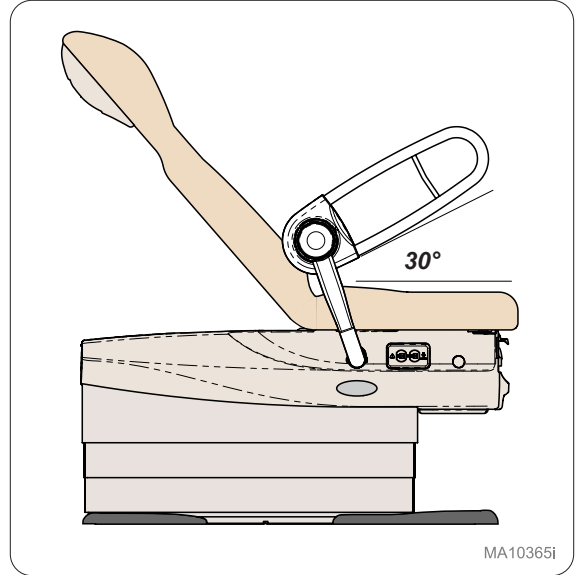
Patient Right - View

The 30° position is used for arm positioning with tall patients.

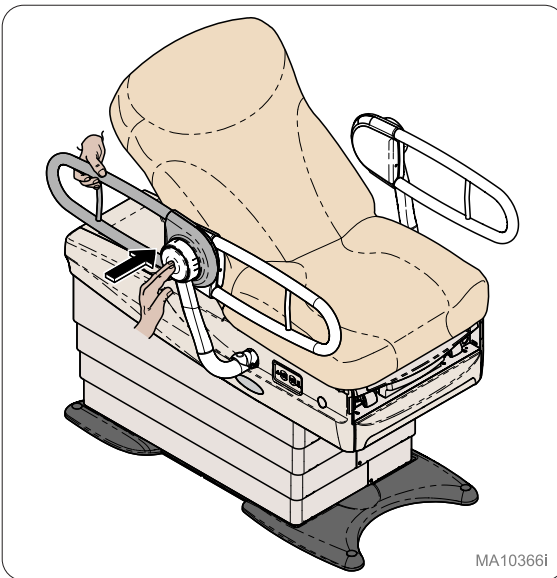Patient Access position


Hold rail with hand while pushing button with other hand to rotate support rail out of way.

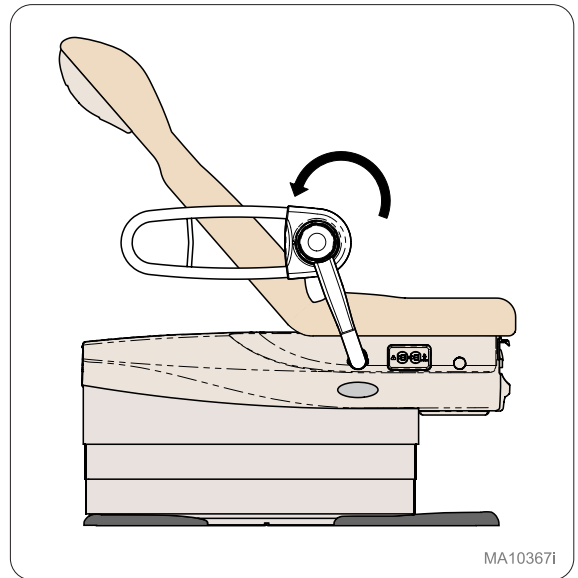
Patient Right - View

This rotation allows for patient access to chair seat.

Operation - Patient Support Rail Plus

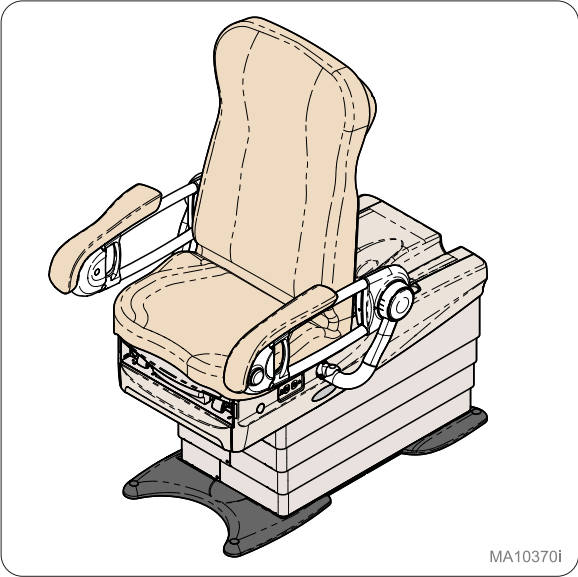
Lowest position
Starting position of Support Rail Plus.

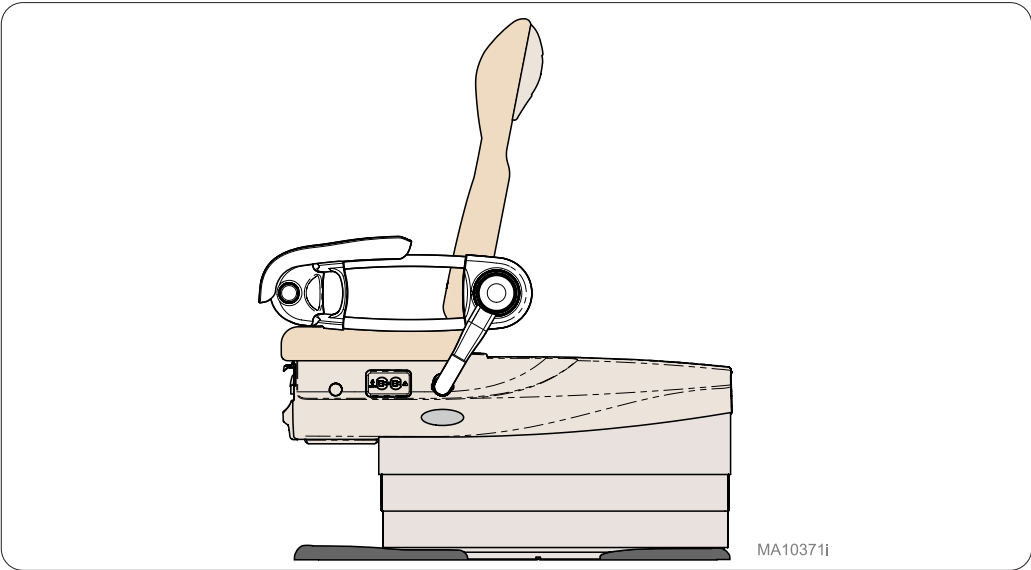Patient Left - View
This position is used for blood pressure checks with short patients.


Operation - Patient Support Rail Plus

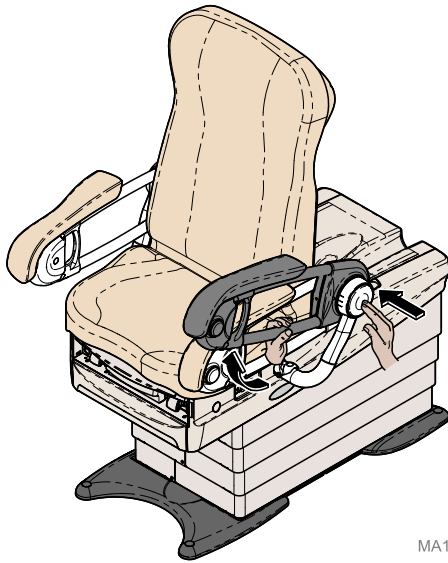MA10372i

Articulation of rail


Hold rail with hand while pushing button with other hand to raise support rail into position.

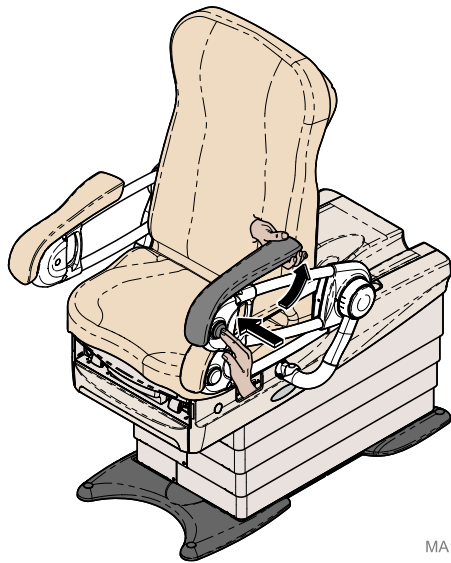MA10373i

Articulation of armrest pad

Hold armrest pad with hand while pushing button with other hand to raise into position.

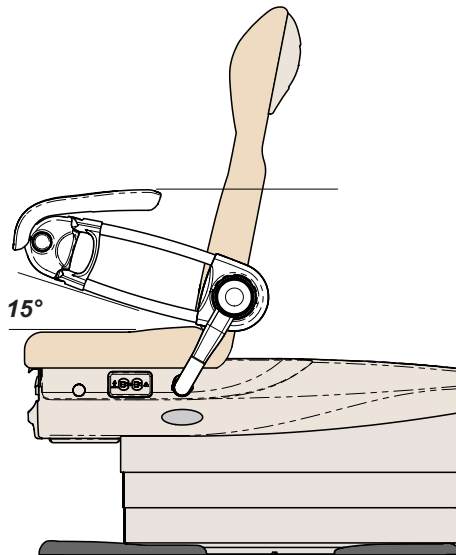
MA10374i

Patient Left - View

The 15° position is used for blood pressure checks with medium height patients.


Operation - Patient Support Rail Plus

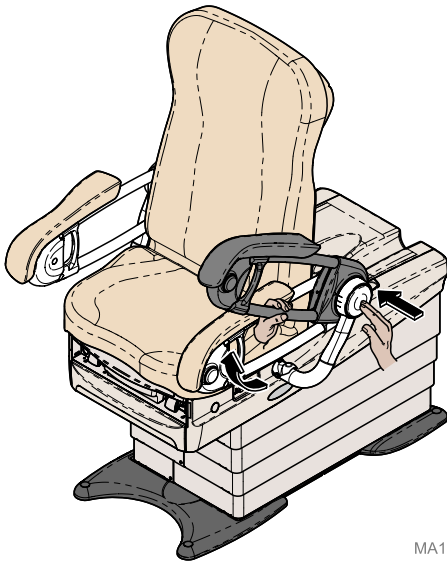MA10375i

Articulation of rail


Hold rail with hand while pushing button with other hand to raise support rail into position.

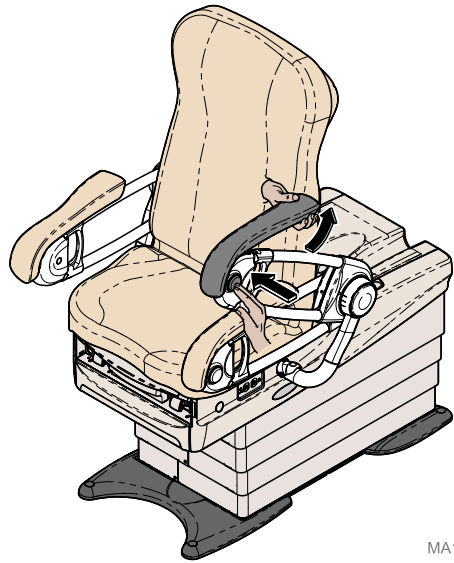MA10376i

Articulation of armrest pad

Hold armrest pad with hand while pushing button with other hand to raise into position.

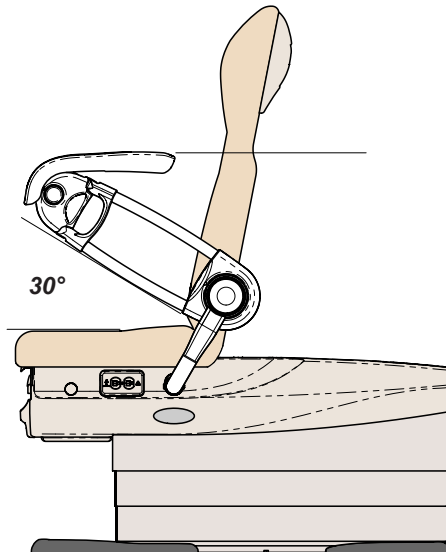
MA10377i

Patient Left - View

The 30° position is used for blood pressure checks with tall patients.

Operation - Patient Support Rail Plus

30°

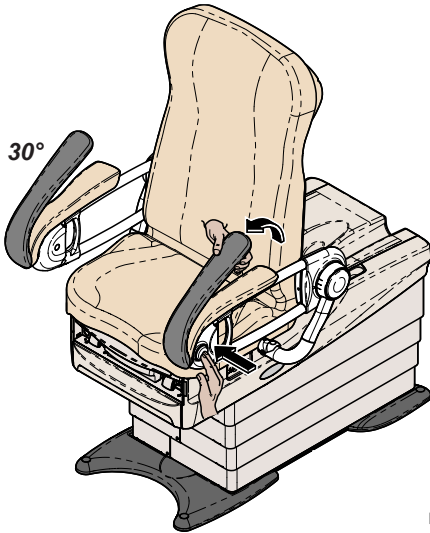
MA10698i

End Transfer position

Hold arm pad with hand while pushing button with other hand to rotate arm pad 30° away from rail.

Ingress Egress Position


30°

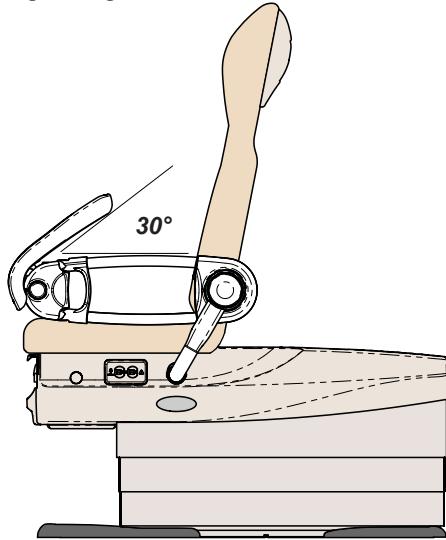MA10699i

Patient Left - View

The ingress egress rotation of 30° for arm pad allows for patients to access chair seat from the front.

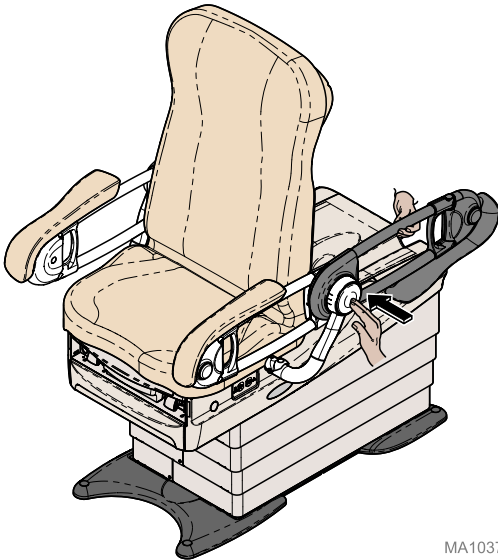MA10378i

Side Transfer position

Hold rail with hand while pushing button with other hand to rotate support rail plus out of way.


Patient Left - View

This rotation allows for patient access to chair seat from the side.


MA10379i

Operation - Patient Support Rail Plus

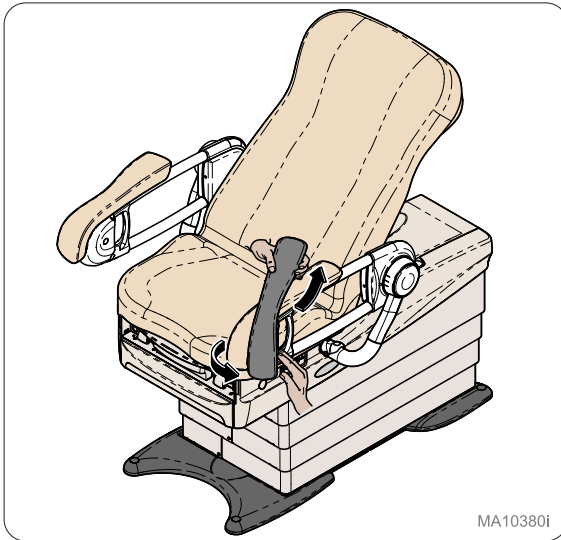
Phlebotomy Articulation

Hold armrest pad with hand while pushing button with other hand to raise into 15° position.

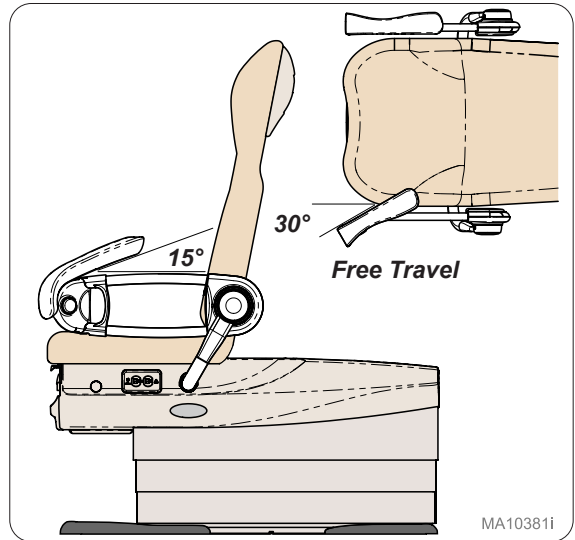
Patient Left - View

This position allows armrest pad to travel side to side freely to position shorter patients arm for phlebotomy procedures.

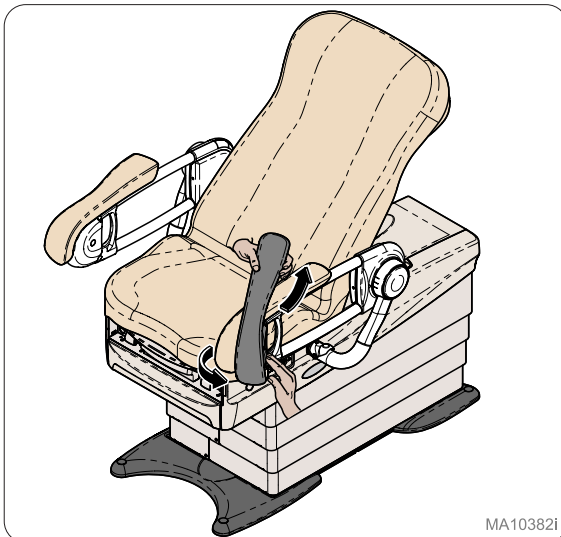Phlebotomy Articulation

Hold armrest pad with hand while pushing button with other hand to raise into 30° position.

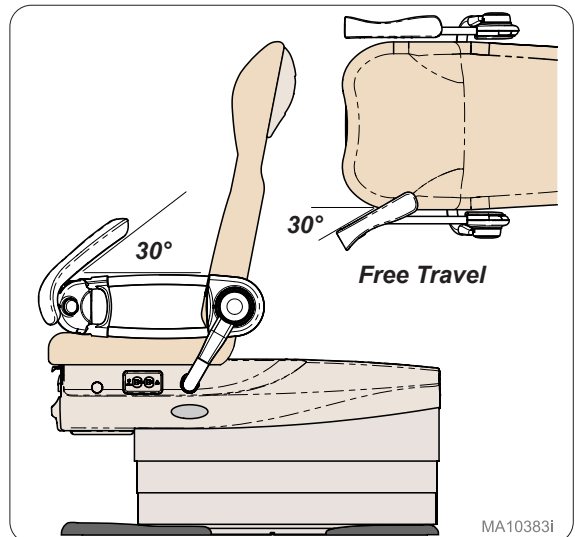
Patient Left - View


This position allows armrest pad to travel side to side freely to position a taller patient's arm for phlebotomy procedures.

Dimensions - Seat Width

Patient Support Rail Plus

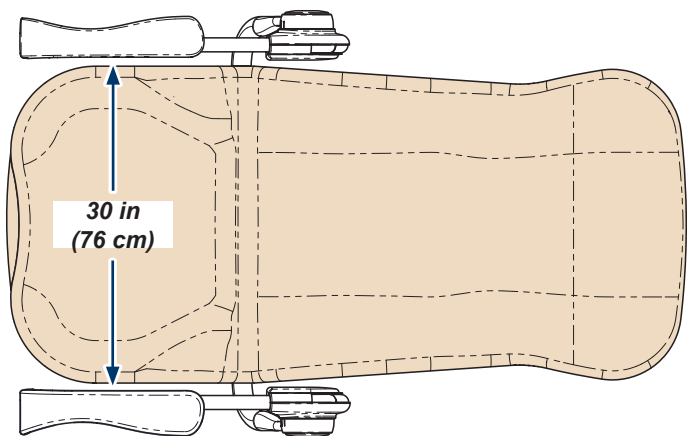Patient Support Rail Plus with Receiver Extension

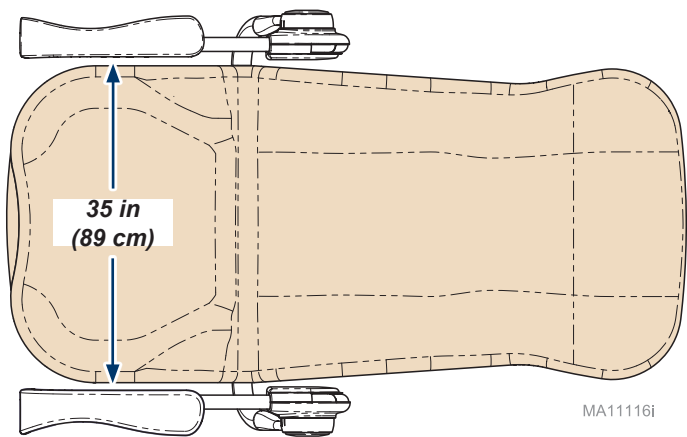MA11116i

Maintenance

Calling for Service

Note

Model / serial number information is required when calling for service.

If service is required, contact Midmark directly:

US 844.856.1230 / Canada 937.526.8585

Monday-Thursday 8:00 AM until 6:00 PM

Friday 8:00 AM until 5:00 PM (ET)

midmark.com

Cleaning

Equipment Alert

The armrest is resistant to most medicinal-type stains, but may be damaged by solvents and dyes. Immediately remove any fluids spilled on the armrest.

Armrest Pad

Wash your armrest weekly with a mild liquid soap and water mixture, rinse with clear water and dry completely to remove disinfection cleaner build-up.

Disinfect your armrest with a solution of standard bleach and water mixed 1 in 10 (10%) or chlorine based cleaners. Follow this with a clear water rinse and thorough drying of material. See current CDC Guideline for Disinfection and Sterilization in Healthcare Facilities.

To minimize damage caused by disinfectant cleaner residue build-up, do not allow disinfectants to pool on the armrest surface. Once the approved contact time has been obtained, remove and dry excess liquid remaining on the surface.

Detailed care and maintenance instructions are included with your product. This information is also available on midmark.com in the Technical library under the User Information tab for your product.

Painted Metal / Plastic Surfaces

Clean the painted metal and plastic surfaces weekly using a clean soft cloth and mild cleaner. Periodic applications of common furniture wax will ease cleaning, and maintain the finished luster of the surfaces.

Periodic Maintenance

Periodically inspect the following areas:

- All fasteners should be in place and tightened securely.
- All mechanical functions should operate properly.

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