


Quick Reference Guide
228 and 229
Barrier-Free® Exam Chair

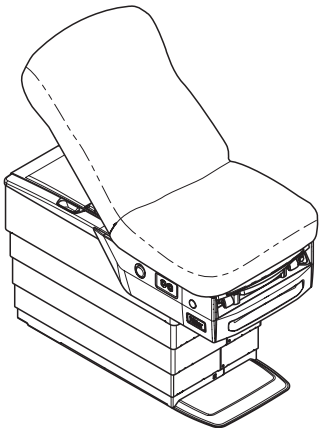

This guide outlines basic operation of the chair. See the User Guide for additional instructions on operation and care, as well as regulatory compliance information.

Model Number:

Serial Number:

Dealer Information:

Purchase Date:

Upholstery Color:

003-10884-00 Rev AA2 (06/01/2026)

Maintenance

Periodically inspect the following areas:

- Power cord must be free of visible damage.
- Foot control cord must be free of visible damage.
- All fasteners should be in place and tightened securely.
- All mechanical functions should operate properly.
- Periodically lubricate the back hinge to maintain quiet, smooth, operation (use light machine oil).
- Have a qualified service technician inspect your equipment every twelve months.

Cleaning

Upholstery/High Touch Metal/Plastic Chair Surfaces

When soiled, follow a two-step process of cleaning and disinfecting as follows. Wash with a mild liquid soap and water mixture, rinse with clean water and completely dry. After it is dry, apply an EPA registered disinfectant. To minimize damage caused by disinfectants, do not allow them to pool on surfaces. Another disinfectant option is to mix standard bleach with water at a ratio of 1:10, respectively. Once the recommended contact time has elapsed, rinse with clean water and dry the disinfected surface.

Alternatively, there are EPA registered products that are both cleaner and disinfectant in one that can be used in a single step per manufacturer instructions.

For additional cleaning and disinfection recommendations, follow the current [Guideline for Disinfection + Sterilization in Healthcare Facilities from the CDC](#).

Detailed care and maintenance instructions are included in the user guide.


Calling for Service

If service is required, contact Midmark directly.

US 844.856.1230 / Canada 937.526.8585
Monday - Thursday 8:00 AM to 6:00 PM ET
Friday 8:00 AM to 5:00 PM ET
[midmark.com](#)

To turn the drawer heater ON / OFF...
Press and release the drawer heater switch.

Note: Switch will illuminate when heater is ON

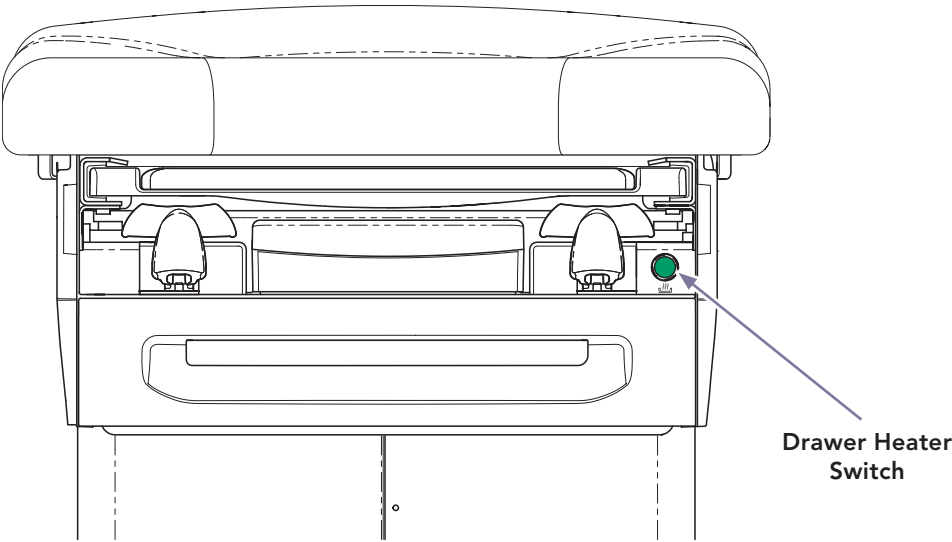

Note

The chair automatically stops at its maximum and minimum height.

To move the chair Up or Down...
Press and hold the appropriate button.



To position the 228 back section...

- A) While holding the the back section in place, squeeze the back release handle (other side).
- B) Move the back section.
- C) Release handle to lock back section in position.

