


Intraoral Dental X-Ray System: Replacing the Horizontal Arm



Caution

The Articulated Arm could spring open and cause serious injury if not secured using the provided tie prior to removal from the horizontal arm. Consult the Installation and Service Manual for guidance on all steps necessary to complete installation.

! CAUTION

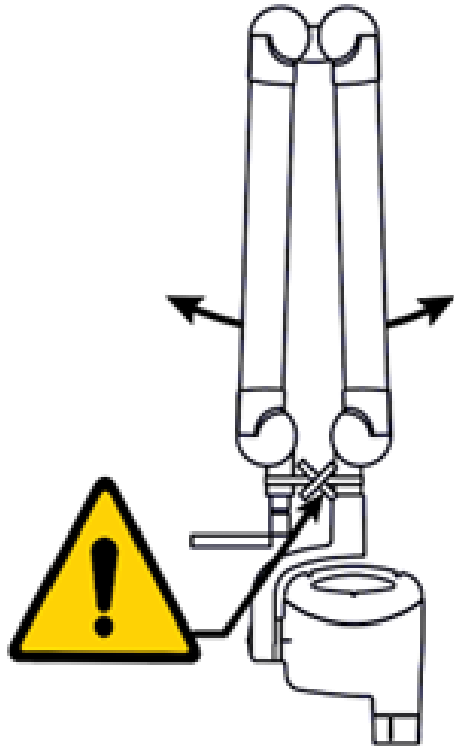

Applies to Models:

DCVx-G2X/H
DCVx-G2X/H/Sx
P701(5 thru 7)
P701(5 thru 8)-P
P701(5 thru 8)V-P
P7017VM
PE701(5 thru 7)
PE701(5 thru 8)-P

Special Tools:

3 mm hex key
10 mm hex socket/wrench
4 mm hex key
7 mm hex socket/wrench
#2 Phillips head screw-
driver

Electronic copies of the Installation and Service Manual can be found at the below URL or QR code.

<https://technicallibrary.midmark.com/WebSite/product-manuals.htm>

