

Paper Roll Holder [9A419]

Components:
Paper Roll Holder ... 1


Special Tools:
none

Intended Use:

To provide a storage location for one roll of table paper.
Max. paper roll size: 18" long x 3.5" dia. (45.7 cm x 9 cm)

Attention:

If necessary, loosen side rail adjustment knobs and pull side rails all the way out. Then tighten knobs.

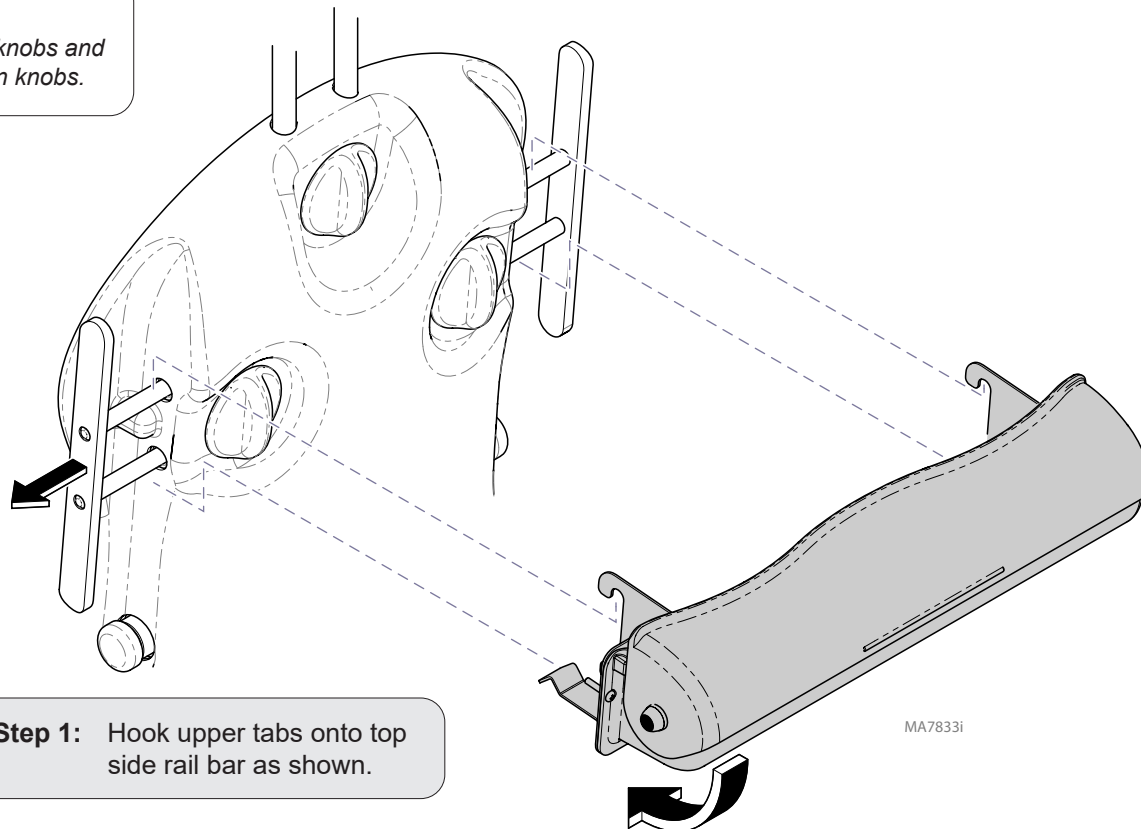

Step 1: Hook upper tabs onto top side rail bar as shown.

Step 2 : Pivot paper roll holder downward until lower tabs snap under bottom side rail bar.